bs-3339R

[Primary Antibody]

BIOSS

www.bioss.com.cn sales@bioss.com.cn techsupport@bioss.com.cn 400-901-9800

phospho-PLCG 2 (Tyr759) Rabbit pAb

- DATASHEET -

Host: Rabbit Isotype: IgG

Clonality: Polyclonal

GenelD: 5336 **SWISS:** P16885

Target: PLCG 2 (Tyr759)

Immunogen: KLH conjugated Synthesised phosphopeptide derived from human

PLC gamma 2 around the phosphorylation site of Tyr759: RM(p-

Y)VD.

Purification: affinity purified by Protein A

Concentration: 1mg/ml

Storage: 0.01M TBS (pH7.4) with 1% BSA, 0.02% Proclin300 and 50%

Glycerol.

Shipped at 4°C. Store at -20°C for one year. Avoid repeated

freeze/thaw cycles.

Background: Enzymes of the phospholipase C family catalyze the hydrolysis of

phospholipids to yield diacylglycerols and water soluble phosphorylated derivatives of the lipid head groups. A number of these enzymes have specificity for phosphoinositides. Of the phosphoinositide specific phospholipase C enzymes, C beta is regulated by heterotrimeric G protein coupled receptors, while the closely related C gamma 1 and C gamma 2 enzymes are controlled by receptor tyrosine kinases. The C gamma 1 and C gamma 2 enzymes are composed of phospholipase domains that flank regions of homology to noncatalytic domains of the SRC oncogene

product, SH2 and SH3.

Applications: IHC-P (1:100-500) **IHC-F** (1:100-500)

IHC-F (1:100-500) **IF** (1:100-500)

Reactivity: Human (predicted: Mouse,


Rat, Rabbit, Pig, Cow, Dog,

Horse)

Predicted MW.: 147 kDa

Subcellular Extracellular matrix ,Cell **Location:** membrane ,Cytoplasm

VALIDATION IMAGES

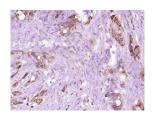

Paraformaldehyde-fixed, paraffin embedded (Human stomach carcinoma); Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15min; Block endogenous peroxidase by 3% hydrogen peroxide for 20 minutes; Blocking buffer (normal goat serum) at 37°C for 30min; Antibody incubation with (Phospho-PLCG 2 (Tyr759)) Polyclonal Antibody, Unconjugated (bs-3339R) at 1:400 overnight at 4°C, followed by operating according to SP Kit(Rabbit) (sp-0023) instructionsand DAB staining.