bs-5811R

[Primary Antibody]


NUAK2 Rabbit pAb



www.bioss.com.cn sales@bioss.com.cn techsupport@bioss.com.cn 400-901-9800

- DATASHEET		400-901-9800	
Host: Rabbit	Isotype: IgG	Applications: IHC-P (1:100-500)	
Clonality: Polyclonal	C C	IHC-F (1:100-500) IF (1:100-500)	
GenelD: 81788	SWISS: Q9H093		
Target: NUAK2		Reactivity: Mouse, Rat (predicted: Human, Rabbit,	
Immunogen: KLH conjugated synthetic peptide derived from human NUAK2: 31-130/628.		Ċow, Dog)	
Purification: affinity purified by Protein A		Predicted MW.: ^{70 kDa}	
Concentration: 1mg/ml		MW.: TO KDa	
Storage: 0.01M TBS (pH7.4) with 1% BSA, 0.02% Proclin300 and 50% Glycerol. Shipped at 4°C. Store at -20°C for one year. Avoid repeated freeze/thaw cycles.		Subcellular Location: ^{Cytoplasm} ,Nucleus	
Background: Stress-activated kinase involved in tolerance to glucose starvation. Induces cell-cell detachment by increasing F-actin conversion to G- actin. Expression is induced by CD95 or TNF-alpha, via NF-kappa-B. Protects cells from CD95-mediated apoptosis and is required for			

– VALIDATION IMAGES

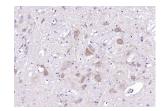

cells.

Paraformaldehyde-fixed, paraffin embedded (mouse cerebellum); Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15min; Block endogenous peroxidase by 3% hydrogen peroxide for 20 minutes; Blocking buffer (normal goat serum) at 37°C for 30min; Antibody incubation with (NUAK2) Polyclonal Antibody, Unconjugated (bs-5811R) at 1:200 overnight at 4°C, followed by operating according to SP Kit(Rabbit) (sp-0023) instructionsand DAB staining.



the increased motility and invasiveness of CD95-activated tumor

Paraformaldehyde-fixed, paraffin embedded (Rat cerebellum); Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15min; Block endogenous peroxidase by 3% hydrogen peroxide for 20 minutes; Blocking buffer (normal goat serum) at 37°C for 30min; Antibody incubation with (NUAK2) Polyclonal Antibody, Unconjugated (bs-5811R) at 1:200 overnight at 4°C, followed by operating according to SP Kit(Rabbit) (sp-0023) instructionsand DAB staining.