bs-21725R

[Primary Antibody]

Bioss

www.bioss.com.cn sales@bioss.com.cn techsupport@bioss.com.cn 400-901-9800

STIM1 Rabbit pAb

- DATASHEET -

Host: Rabbit Isotype: IgG

Clonality: Polyclonal

GeneID: 6786 **SWISS:** Q13586

Target: STIM1

Immunogen: KLH conjugated synthetic peptide derived from human STIM1:

601-685/685.

Purification: affinity purified by Protein A

Concentration: 1mg/ml

Storage: 0.01M TBS (pH7.4) with 1% BSA, 0.02% Proclin300 and 50%

Glycerol.

Shipped at 4°C. Store at -20°C for one year. Avoid repeated

freeze/thaw cycles.

Background: Plays a role in mediating Ca(2+) influx following depletion of

intracellular Ca(2+) stores. Acts as Ca(2+) sensor in the endoplasmic reticulum via its EF-hand domain. Upon Ca(2+) depletion, translocates from the endoplasmic reticulum to the

plasma membrane where it activates the Ca(2+) release-activated

Ca(2+) (CRAC) channel subunit, TMEM142A/ORAI1.

Applications: IHC-P (1:100-500)

IHC-F (1:100-500) **IF** (1:100-500)

Reactivity: Human, Rat


(predicted: Mouse, Sheep,

Cow, Horse)


Predicted MW.: 75 kDa

Subcellular Cell membrane Cytoplasm

VALIDATION IMAGES

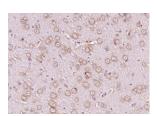

Paraformaldehyde-fixed, paraffin embedded (rat brain); Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15min; Block endogenous peroxidase by 3% hydrogen peroxide for 20 minutes; Blocking buffer (normal goat serum) at 37°C for 30min; Antibody incubation with (STIM1) Polyclonal Antibody, Unconjugated (bs-21725R) at 1:200 overnight at 4°C, followed by operating according to SP Kit(Rabbit) (sp-0023) instructions and DAB staining.

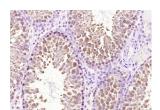

Paraformaldehyde-fixed, paraffin embedded (rat testis); Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15min; Block endogenous peroxidase by 3% hydrogen peroxide for 20 minutes; Blocking buffer (normal goat serum) at 37°C for 30min; Antibody incubation with (STIM1) Polyclonal Antibody, Unconjugated (bs-21725R) at 1:200 overnight at 4°C, followed by operating according to SP Kit(Rabbit) (sp-0023) instructionsand DAB staining.