## bs-8574R

# [ Primary Antibody ]

# Bioss

www.bioss.com.cn sales@bioss.com.cn techsupport@bioss.com.cn 400-901-9800

# ox-LDL Rabbit pAb

- DATASHEET -

Host: Rabbit Isotype: IgG

Clonality: Polyclonal

Target: ox-LDL

**Purification:** affinity purified by Protein A

Concentration: 1mg/ml

Storage: 0.01M TBS (pH7.4) with 1% BSA, 0.02% Proclin300 and 50%

Glycerol.

Shipped at 4°C. Store at -20°C for one year. Avoid repeated

freeze/thaw cycles.

**Background:** Low-density lipoprotein (LDL) is the carrier protein for cholesterol

in the blood. LDL binds to its receptor on the capillary walls and thereby mediates the uptake and clearence of cholesterol from the circulation. In atherosclerotic lesions oxidatively modified LDL is found and oxidized LDL is specifically recognized and ingested by macrophages via scavenger receptor A and CD36. Oxidized LDL may be a marker of atherosclerosis but the precise changes in oxidized LDL are not well described. Low-density lipoprotein

oxidised with Cu2SO4.

Applications: IHC-P (1:100-500)

IHC-F (1:100-500) IF (1:100-500) ICC/IF (1:100)

Reactivity: Human (predicted: Mouse)


Predicted MW.: 31 kDa

**Subcellular Location:** Secreted ,Cell membrane


### VALIDATION IMAGES -



Paraformaldehyde-fixed, paraffin embedded Human Liver; Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15 min; Antibody incubation with ox-LDL Polyclonal Antibody, Unconjugated (bs-8574R) at 1:200 overnight at 4°C, followed by conjugation to the bs-0295G-HRP and DAB (C-0010) staining.

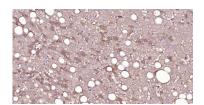

Paraformaldehyde-fixed, paraffin embedded Human Liver Cancer; Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15 min; Antibody incubation with ox-LDL Polyclonal Antibody, Unconjugated (bs-8574R) at 1:200 overnight at 4°C, followed by conjugation to the bs-0295G-HRP and DAB (C-0010) staining.

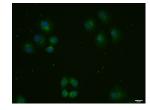

HepG2 cell; 4% Paraformaldehyde-fixed; Triton X-100 at room temperature for 20 min; Blocking buffer (normal goat serum, C-0005) at 37°C for 20 min; Antibody incubation with (ox-LDL) polyclonal Antibody, Unconjugated (bs-8574R) 1:100, 90 minutes at 37°C; followed by a conjugated Goat Anti-Rabbit IgG antibody at 37°C for 90 minutes, DAPI (blue, C02-04002) was used to stain the cell nuclei.

### - SELECTED CITATIONS -

- [IF=3.23] Jehle J, Hoyer FF, Schöne B, Pfeifer P, Schild K, Jenniches I, et al. (2016) Myeloid-Specific Deletion of Diacylglycerol Lipase α Inhibits Atherogenesis in ApoE-Deficient Mice. PLoS ONE 11(1): e0146267. IHC;="Mouse". 26731274
- [IF=1.58] Liu, Liang, et al. "Quercetin Alleviates High-Fat Diet-Induced Oxidized Low-Density Lipoprotein Accumulation in the Liver: Implication for Autophagy Regulation." BioMed Research International 2015 (2015). IHC;="Mouse". 26697490