bs-0011R

[Primary Antibody]

www.bioss.com.cn sales@bioss.com.cn techsupport@bioss.com.cn 400-901-9800

Gastrin Releasing Peptide Rabbit pAb

DATASHEET -

Host: Rabbit Isotype: IgG

Clonality: Polyclonal

GeneID: 2922 SWISS: P07492

Target: Gastrin Releasing Peptide

Immunogen: KLH conjugated synthetic peptide derived from human GRP:

Purification: affinity purified by Protein A

Concentration: 1mg/ml

Storage: 0.01M TBS (pH7.4) with 1% BSA, 0.02% Proclin300 and 50%

Glycerol.

Shipped at 4°C. Store at -20°C for one year. Avoid repeated

freeze/thaw cycles.

Background: This gene encodes a member of the bombesin-like family of gastrin-releasing peptides. The encoded preproprotein is proteolytically processed to generate two peptides, gastrinreleasing peptide and neuromedin-C. These peptides regulate numerous functions of the gastrointestinal and central nervous systems, including release of gastrointestinal hormones, smooth muscle cell contraction, and epithelial cell proliferation. These peptides are also likely to play a role in human cancers of the lung, colon, stomach, pancreas, breast, and prostate. Alternative splicing results in multiple transcript variants, at least one of which encodes a preproprotein that is proteolytically processed. [provided by RefSeq, Jan 2016].

Applications: IHC-P (1:100-500)


IHC-F (1:100-500) **IF** (1:100-500)

Reactivity: Human

Predicted 3/14 kDa MW.:

Subcellular Location: Secreted

VALIDATION IMAGES

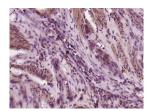

Paraformaldehyde-fixed, paraffin embedded (human gastric carcinoma); Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15min; Block endogenous peroxidase by 3% hydrogen peroxide for 20 minutes; Blocking buffer (normal goat serum) at 37°C for 30min; Antibody incubation with (Gastrin Releasing Peptide) Polyclonal Antibody, Unconjugated (bs-0011R) at 1:400 overnight at 4°C, followed by operating according to SP Kit(Rabbit) (sp-0023) instructions and DAB staining.



Paraformaldehyde-fixed, paraffin embedded (human brain glioma); Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15min; Block endogenous peroxidase by 3% hydrogen peroxide for 20 minutes; Blocking buffer (normal goat serum) at 37°C for 30min; Antibody incubation with (Gastrin Releasing Peptide) Polyclonal Antibody, Unconjugated (bs-0011R) at 1:400 overnight at 4°C, followed by operating according to SP Kit(Rabbit) (sp-0023) instructions and DAB staining.