

**bs-16516R****[ Primary Antibody ]****ZNF676 Rabbit pAb****Bioss**  
**ANTIBODIES**

www.bioss.com.cn


sales@bioss.com.cn

techsupport@bioss.com.cn


400-901-9800

**— DATASHEET —**


<b>Host:</b> Rabbit	<b>Isotype:</b> IgG	<b>Applications:</b> <b>WB</b> (1:500-2000) <b>IHC-P</b> (1:100-500) <b>IHC-F</b> (1:100-500) <b>IF</b> (1:100-500)  <b>Reactivity:</b> Human  <b>Predicted MW.:</b> 68 kDa  <b>Subcellular Location:</b> Nucleus
<b>Clonality:</b> Polyclonal		
<b>GeneID:</b> 163223	<b>SWISS:</b> Q8N7Q3	
<b>Target:</b> ZNF676		
<b>Immunogen:</b> KLH conjugated synthetic peptide derived from human ZNF676: 501-588/588.		
<b>Purification:</b> affinity purified by Protein A		
<b>Concentration:</b> 1mg/ml		
<b>Storage:</b> 0.01M TBS (pH7.4) with 1% BSA, 0.02% Proclin300 and 50% Glycerol. Shipped at 4°C. Store at -20°C for one year. Avoid repeated freeze/thaw cycles.		

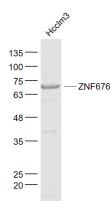**— VALIDATION IMAGES —**

Sample: Hcclm3(Human) Cell Lysate at 30 ug  
Primary: Anti- ZNF676 (bs-16516R) at 1/300  
dilution Secondary: IRDye800CW Goat Anti-  
Rabbit IgG at 1/20000 dilution Predicted band  
size: 68 kD Observed band size: 68 kD

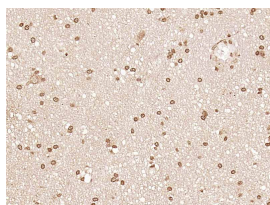

Paraformaldehyde-fixed, paraffin embedded  
(Human brain glioma); Antigen retrieval by  
boiling in sodium citrate buffer (pH6.0) for  
15min; Block endogenous peroxidase by 3%  
hydrogen peroxide for 20 minutes; Blocking  
buffer (normal goat serum) at 37°C for 30min;  
Antibody incubation with (ZNF676) Polyclonal  
Antibody, Unconjugated (bs-16516R) at 1:400  
overnight at 4°C, followed by operating  
according to SP Kit(Rabbit) (sp-0023)  
instructions and DAB staining.

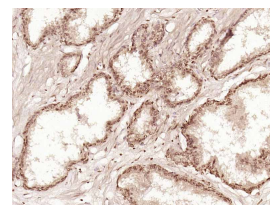

Paraformaldehyde-fixed, paraffin embedded  
(Human prostate); Antigen retrieval by boiling in  
sodium citrate buffer (pH6.0) for 15min; Block  
endogenous peroxidase by 3% hydrogen  
peroxide for 20 minutes; Blocking buffer (normal  
goat serum) at 37°C for 30min; Antibody  
incubation with (ZNF676) Polyclonal Antibody,  
Unconjugated (bs-16516R) at 1:400 overnight at  
4°C, followed by operating according to SP  
Kit(Rabbit) (sp-0023) instructions and DAB  
staining.