

## CD68 Ready-To-Use IHC Kit

Cat.No: IHC0133  
Applications: **IHC-P**  
Reactivity: Human, Rat  
Size: 50T  
Assay type: Immunohistochemistry  
Sample type: FFPE tissue  
General Information:

Number	Component	Size	Concentration	Storage
1	PBS Buffer (powder)	2 L X 2	20x	RT
2	Antigen Retrieval Buffer	20 ml	100x	2-8°C
3	Endogenous Peroxidase Blocking Buffer	3 ml	RTU	2-8°C, protect from light
4	Blocking Buffer	3 ml	RTU	2-8°C
5	Primary Antibody (CD68 Mouse mAb)	6 ml	RTU	2-8°C
6	Secondary Antibody (HRP-Goat anti-Mouse IgG pAb)	6 ml	RTU	2-8°C
7	Chromogen Component A	0.3 ml	RTU	-20°C, protect from light
8	Chromogen Component B	0.3 ml	RTU	-20°C
9	Counter Staining Reagent	5 ml	RTU	RT
10	Mounting Media	5 ml	RTU	RT
11	Control slide (大鼠脾, 人脾)	2 slides	RTU	RT
12	Datasheet	1 copy		

### Immunohistochemistry Protocol:

#### 1. Deparaffinization And Rehydration

Immerse slides in fresh xylene for 15 minutes and then repeat two more times using separate containers. Immerse slides sequentially in 100%, 95%, 90%, 80%, and 70% ethanol solutions for 5 minutes each. Rinse slides 3 times with distilled water for 5 minutes each.

#### 2. Antigen Retrieval (Pressure Cooker)

Prepare a 1x antigen retrieval solution by diluting the 100x Antigen Retrieval Buffer using distilled water. Add the appropriate amount of 1x antigen retrieval solution into the pressure cooker and place a heat-resistant staining container filled with the same solution inside the cooker. Heat the solution to boiling with the lid of the pressure cooker rested on top without being secured. Once it's boiling, transfer the slides from the distilled water to the staining container inside the pressure cooker. Follow the manufacturer's instructions to secure the lid of the pressure cooker. As soon as the cooker reaches full pressure, time

three minutes. After three minutes, move the pressure cooker to an empty sink and cool it down by running cold water over the cooker. Once depressurized, open the lid and transfer the staining container with the slides to room temperature. After 20 minutes, rinse 3 times with PBS Buffer (dissolve the powder in 2L distilled water) for 5 minutes each.

### 3. **Block Endogenous Peroxidase**

Drain the liquid off the slides and then use a hydrophobic IHC pen to draw circles on the slides around tissue sections. Add 2-4 drops of **Endogenous Peroxidase Blocking Buffer** directly on slides, covering the whole tissue and block slides for 15 minutes at RT. Rinse 3 times with **PBS Buffer** for 5 minutes each.

### 4. **Serum Blocking**

Block with 2-4 drops of **Blocking Buffer** for 20 minutes at RT.

### 5. **Primary Antibody Incubation**

Drain blocking buffer from slides. Incubate slides with 2-4 drops of **CD68 Mouse mAb** overnight at 4°C or 1-2 hours at RT. Rinse 3 times with **PBS Buffer** for 5 minutes each.

### 6. **Secondary Antibody Incubation**

Incubate slides with 2-4 drops of **HRP-Goat anti-Mouse IgG pAb** for 1-2 hours at RT. Rinse slides 3 times with **PBS Buffer** for 5 minutes each.

### 7. **Signal Development**

Remove residual liquid around the tissue section. Add 50ul fresh **DAB Buffer (Chromogen Component A : Chromogen Component B : PBS Buffer=1:1:18)** to cover the tissue. Monitor the reaction under the microscope until a brown color is visible (approximate 3-5 minutes at RT). Stop reaction immediately by rinsing with distilled water. Rinse slides 3 times with distilled water for 5 minutes each.

### 8. **Counterstain**

Counterstain with an appropriate amount of **Counter Staining Reagent** for 3-5 minutes at RT. Rinse slides with distilled water for 5 minutes. Use 2-4 drops of **Differentiation reagent** to cover the tissue for 30 seconds. Rinse slides twice with distilled water for 5 minutes each.

### 9. **Dehydration Sheet**

Immerse slides sequentially in 70%, 80%, 90%, 95%, and 100% ethanol for 5 minutes each at RT. Immerse slides in 2 changes of fresh xylene, 15 minutes each. Drop some **Mounting Media** on the tissue. Mount coverslips.

## Notes:


1. The positive control slide provided in the kit allows you to be sure that the experimental set-up is working properly.
2. Do not allow slides to dry at any time during this procedure.
3. Please don't replace the matching reagents in this product with other manufacturers' products.
4. As DAB is a carcinogen, please take necessary precautions.
5. PBS (reagent 1) can be stored for one week at 4°C after preparation; The antigen retrieval buffer (1× reagent 2) and the chromogenic agent (the mixture of reagents 7 and 8) should be prepared right before each assay.

***Please cite this product as "IHC0133, Bioss Antibodies". Citation example: "Tissue sections using CD68 IHC Kit (IHC0133, Bioss Antibodies) were stained for CD68 according to the manufacturer's instructions."***


## Introduction:

CD68 (Macrosialin) is a 110 kDa integral membrane glycoprotein predominantly expressed on the intracellular lysosomes of monocytes and macrophages and to a lesser extent by dendritic cells and peripheral blood granulocytes. Also, CD68 could play a role in phagocytic activities of tissue macrophages, both in intracellular lysosomal metabolism and extracellular cell-cell and cell-pathogen interactions. CD68 is expressed by interdigitating reticulum cells in tonsil and some histiocytic lymphoma or histiocytosis, acute myeloid leukemia (AML), and granulocytic sarcoma. Elevated expression of CD68 has been demonstrated on CD34+ cells in various human malignancies, including several Acute Myeloid Leukemia studies.


## Validation Data

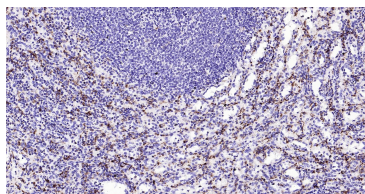


Immunohistochemical analysis of paraffin embedded human spleen tissue slide using IHC0133 (CD68 IHC Kit).

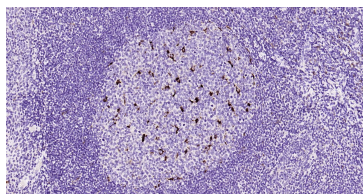

Immunohistochemical analysis of paraffin embedded human tonsil tissue slide using IHC0133 (CD68 IHC Kit).

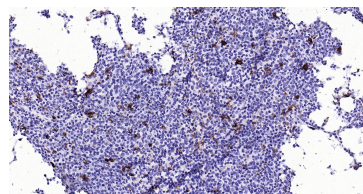

Immunohistochemical analysis of paraffin embedded human diffuse large B-cell lymphoma tissue slide using IHC0133 (CD68 IHC Kit).

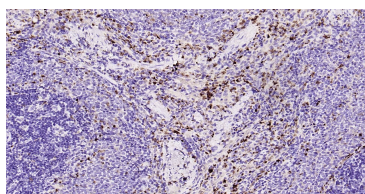

Immunohistochemical analysis of paraffin  
embedded rat spleen tissue slide using IHC0133  
(CD68 IHC Kit).