

www.bioss.com.cn sales@bioss.com.cn techsupport@bioss.com.cn 400-901-9800

Rat Lamin B Ready-To-Use IHC Kit

Cat.No: IHC0183R

Applications: IHC-P
Reactivity: Rat
Size: 50T

Assay type: Immunohistochemistry

Sample type: FFPE tissue

General Information:

Number	Component	Size	Concentration	Storage
1	PBS Buffer (powder)	2L×2	20x	RT
2	Antigen Retrieval Buffer	20 ml	100x	2-8°C
3	Endogenous Peroxidase Blocking Buffer	3 ml	RTU	2-8°C, protect from light
4	Blocking Buffer	3 ml	RTU	2-8°C
5	Primary Antibody (Rat Lamin B (Nuclear Loading Control) Mouse mAb)	6 ml	RTU	2-8°C
6	Secondary Antibody (HRP-Goat anti-Mouse IgG pAb)	6 ml	RTU	2-8°C
7	Chromogen Component A	0.3 ml	RTU	-20°C,protect from light
8	Chromogen Component B	0.3 ml	RTU	-20°C
9	Counter Staining Reagent	5 ml	RTU	RT
10	Mounting Media	5 ml	RTU	RT
11	Control slide (Rat brain)	1 slide	RTU	RT
12	Datasheet	1 сору		

Immunohistoche mistry Protocol:

1. Deparaffinization And Rehydration

Immerse slides in fresh xylene for 15 minutes and then repeat two more times using separate containers. Immerse slides sequentially in 100%, 95%, 90%, 80%, and 70% ethanol solutions for 5 minutes each. Rinse slides 3 times with distilled water for 5 minutes each.

2. Antigen Retrieval

Add $100 \times$ Antigen Retrieval Buffer into distilled water to prepare a $1 \times$ solution. Boil slides in $1 \times$ solution at 95° C- 100° C for 15 minutes. Move the slides to $1 \times$ solution at room temperature (RT) and allow them to stand for 20 minutes. Rinse 3 times with **PBS Buffer** (dissolve the powder in 2L distilled water) for 5 minutes each.

3. Block Endogenous Peroxidase

Drain the liquid off the slides and then use a hydrophobic IHC pen to draw circles on the slides around tissue sections. Add 2-4 drops of **Endogenous Peroxidase Blocking Buffer** directly on slides, covering the whole tissue and block slides for 15 minutes at RT. Rinse 3 times with **PBS Buffer** for 5 minutes each.

4. Serum Blocking

Block with 2-4 drops of **Blocking Buffer** for 20 minutes at RT.

5. Primary Antibody Incubation

Drain blocking buffer from slides. Incubate slides with 2-4 drops of Rat Lamin B (Nuclear Loading Control) Mouse mAb overnight at 4°C or 1-2 hours at RT. Rinse 3 times with PBS Buffer for 5 minutes each.

6. Secondary Antibody Incubation

Incubate slides with 2-4 drops of **HRP-Goat anti-Mouse IgG pAb** for 1-2 hours at RT. Rinse slides 3 times with **PBS Buffer** for 5 minutes each.

7. Signal Development

Remove residual liquid around the tissue section. Add 50ul fresh **DAB Buffer** (**Chromogen Component A : Chromogen Component B : PBS Buffer=1:1:18**) to cover the tissue. Monitor the reaction under the microscope until a brown color is visible (approximate 3-5 minutes at RT). Stop reaction immediately by rinsing with distilled water. Rinse slides 3 times with distilled water for 5 minutes each.

8. Counterstain

Counterstain with an appropriate amount of **Counter Staining Reagent** for 3-5 minutes at RT. Rinse slides with distilled water for 5 minutes. Use 2-4 drops of **Differentiation reagent** to cover the tissue for 30 seconds. Rinse slides twice with distilled water for 5 minutes each.

9. Dehydration Sheet

Immerse slides sequentially in 70%, 80%, 90%, 95%, and 100% ethanol for 5 minutes each at RT. Immerse slides in 2 changes of fresh xylene, 15 minutes each. Drop some **Mounting**Media on the tissue. Mount coverslips.

Notes:

- 1. The positive control slide provided in the kit allows you to be sure that the experimental set-up is working properly.
- 2. Do not allow slides to dry at any time during this procedure.
- 3. Please don't replace the matching reagents in this product with other

manufacturers' products.

4. As DAB is a carcinogen, please take necessary precautions.


5. PBS (reagent 1) can be stored for one week at 4°C after preparation; The antigen retrieval buffer (1×reagent 2) and the chromogenic agent (the mixture of reagents 7 and 8) should be prepared right before each assay.

Please cite this product as "IHC0183R, Bioss Antibodies". Citation example: "Rat Tissue sections using LMNB1 IHC Kit (IHC0183R, Bioss Antibodies) were stained for LMNB1 according to the manufacturer's instructions."

Introduction:

Laminins, a family of extracellular matrix glycoproteins, are the major noncollagenous constituent of basement membranes. They have been implicated in a wide variety of biological processes including cell adhesion, differentiation, migration, signaling, neurite outgrowth and metastasis. Laminins are composed of 3 non identical chains: laminin alpha, beta and gamma (formerly A, B1, and B2, respectively) and they form a cruciform structure consisting of 3 short arms, each formed by a different chain, and a long arm composed of all 3 chains. Defects in the LMNB1 gene can cause Leukodystrophy, demyelinating, autosomal dominant, adult-onset (ADLD).

Validation Data

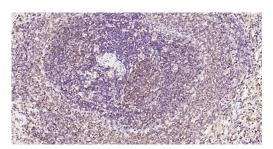

Immunohistochemical analysis of paraffin embedded rat spleen tissue slide using IHC0183R (Rat Lamin B IHC Kit).