bs-16929R

- DATASHEET -

[Primary Antibody]


KCTD17 Rabbit pAb




www.bioss.com.cn sales@bioss.com.cn techsupport@bioss.com.cn 400-901-9800

- DATASHEET		
Host: Rabbit	lsotype: lgG	Applications: IHC-P (1:100-500)
Clonality: Polyclonal		IHC-F (1:100-500) IF (1:100-500)
GenelD: 79734	SWISS: Q8N5Z5	(,
Target: KCTD17		Reactivity: Mouse, Rat (predicted: Human, Pig,
Immunogen: KLH conjugated synthetic peptide derived from human KCTD17: 1-100/521.		Cow)
Purification: affinity purified by Protein A		Predicted MW.: ^{36 kDa}
Concentration: 1mg/ml		MW.: CONSU
Storage: 0.01M TBS (pH7.4) with 1% BSA, 0.02% Proclin300 and 50% Glycerol. Shipped at 4°C. Store at -20°C for one year. Avoid repeated freeze/thaw cycles.		Subcellular Location: Cytoplasm
Background: KCTD16 (potassium channel tetramerisation domain containing 16), also known as BTB/POZ domain-containing protein KCTD16, is a 428 amino acid protein that contains one BTB (POZ) domain. An auxiliary subunit of GABAB R1 and GABAB R2, KCTD16 increases		

- VALIDATION IMAGES

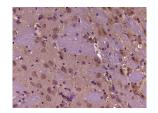

Paraformaldehyde-fixed, paraffin embedded (Rat brain); Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15min; Block endogenous peroxidase by 3% hydrogen peroxide for 20 minutes; Blocking buffer (normal goat serum) at 37°C for 30min; Antibody incubation with (KCTD17) Polyclonal Antibody, Unconjugated (bs-16929R) at 1:400 overnight at 4°C, followed by operating according to SP Kit(Rabbit) (sp-0023) instructions and DAB staining.

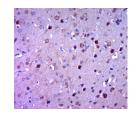

agonist potency and alters the G-protein signaling of the receptors

by accelerating onset and promoting desensitization.

Paraformaldehyde-fixed, paraffin embedded (mouse brain tissue); Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15min; Block endogenous peroxidase by 3% hydrogen peroxide for 20 minutes; Blocking buffer (normal goat serum) at 37°C for 30min; Antibody incubation with (KCTD17) Polyclonal Antibody, Unconjugated (bs-16929R) at 1:400 overnight at 4°C, followed by a conjugated secondary (sp-0023) for 20 minutes and DAB staining.

Important Note: This product as supplied is intended for research use only, not for use in human, therapeutic or diagnostic applications.