

www.bioss.com.cn sales@bioss.com.cn techsupport@bioss.com.cn 400-901-9800

Ran Rabbit pAb

Catalog Number: bs-4282R

Target Protein: Ran

Concentration: 1mg/ml

Form: Liquid

Host: Rabbit

Clonality: Polyclonal

Isotype: IgG

Applications: IHC-P (1:100-500), IHC-F (1:100-500), IF (1:100-500)

Reactivity: Human, Rat (predicted:Mouse, Pig, Cow, Chicken, Dog)

Predicted MW: 26 kDa Entrez Gene: 5901 Swiss Prot: P62826

Source: KLH conjugated synthetic peptide derived from human Ran: 120-216/216.

Purification: affinity purified by Protein A

Storage: 0.01M TBS (pH7.4) with 1% BSA, 0.02% Proclin300 and 50% Glycerol.


Shipped at 4°C. Store at -20°C for one year. Avoid repeated freeze/thaw cycles.

Background: GTP-binding protein involved in nucleocytoplasmic transport. Required for the import of


protein into the nucleus and also for RNA export. Involved in chromatin condensation and control of cell cycle (By similarity). The complex with BIRC5/ survivin plays a role in mitotic spindle formation by serving as a physical scaffold to help deliver the RAN effector molecule

TPX2 to microtubules.

VALIDATION IMAGES

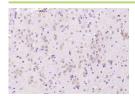

Paraformaldehyde-fixed, paraffin embedded (rat brain); Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15min; Block endogenous peroxidase by 3% hydrogen peroxide for 20 minutes; Blocking buffer (normal goat serum) at 37°C for 30min; Antibody incubation with (Ran) Polyclonal Antibody, Unconjugated (bs-4282R) at 1:200 overnight at 4°C, followed by operating according to SP Kit(Rabbit) (sp-0023) instructions and DAB staining.

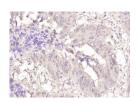

Paraformaldehyde-fixed, paraffin embedded (human rectal carcinoma); Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15min; Block endogenous peroxidase by 3% hydrogen peroxide for 20 minutes; Blocking buffer (normal goat serum) at 37°C for 30min; Antibody incubation with (Ran) Polyclonal Antibody, Unconjugated (bs-4282R) at 1:200 overnight at 4°C, followed by operating according to SP Kit(Rabbit) (sp-0023) instructions and DAB staining.