bs-17116R

[Primary Antibody]

www.bioss.com.cn sales@bioss.com.cn techsupport@bioss.com.cn 400-901-9800

RNF109/TRIM56 Rabbit pAb

DATASHEET -

Host: Rabbit Isotype: IgG

Clonality: Polyclonal

GenelD: 81844 SWISS: Q9BRZ2

Target: RNF109/TRIM56

Immunogen: KLH conjugated synthetic peptide derived from human

RNF109/TRIM56: 651-755/755.

Purification: affinity purified by Protein A

Concentration: 1mg/ml

Storage: 0.01M TBS (pH7.4) with 1% BSA, 0.02% Proclin300 and 50%

Glycerol.

Shipped at 4°C. Store at -20°C for one year. Avoid repeated

freeze/thaw cycles.

Background: E3 ubiquitin-protein ligase that mediates 'Lys-63'-linked

polyubiquitination of TMEM173/STING, thereby playing a key role in innate immunity. TMEM173/STING 'Lys-63'-linked ubiquitination activates the production of type I interferon IFN-beta following detection of pathogen- and host-derived double-stranded DNA.

Applications: IHC-P (1:100-500)

IHC-F (1:100-500) **IF** (1:100-500)

Reactivity: Rat (predicted: Human,


Mouse, Rabbit, Sheep, Cow,

Dog, Horse)

Predicted 81 kDa MW.:

Subcellular Location: Secreted

VALIDATION IMAGES

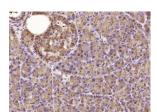

Paraformaldehyde-fixed, paraffin embedded (rat pancreas); Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15min; Block endogenous peroxidase by 3% hydrogen peroxide for 20 minutes; Blocking buffer (normal goat serum) at 37°C for 30min; Antibody incubation with (RNF109 TRIM56) Polyclonal Antibody, Unconjugated (bs-17116R) at 1:200 overnight at 4°C, followed by operating according to SP Kit(Rabbit) (sp-0023) instructions and DAB staining.



Paraformaldehyde-fixed, paraffin embedded (rat thyroid gland); Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15min; Block endogenous peroxidase by 3% hydrogen peroxide for 20 minutes; Blocking buffer (normal goat serum) at 37°C for 30min; Antibody incubation with (RNF109 TRIM56) Polyclonal Antibody, Unconjugated (bs-17116R) at 1:200 overnight at 4°C, followed by operating according to SP Kit(Rabbit) (sp-0023) instructions and DAB staining.