

www.bioss.com.cn sales@bioss.com.cn techsupport@bioss.com.cn 400-901-9800

## RNF109/TRIM56 Rabbit pAb

Catalog Number: bs-17116R

Target Protein: RNF109/TRIM56

Concentration: 1mg/ml

Form: Liquid Host: Rabbit

Clonality: Polyclonal

Isotype: IgG

Applications: IHC-P (1:100-500), IHC-F (1:100-500), IF (1:100-500)

Reactivity: Rat (predicted: Human, Mouse, Rabbit, Sheep, Cow, Dog, Horse)

Predicted MW: 81 kDa
Entrez Gene: 81844
Swiss Prot: Q9BRZ2

Source: KLH conjugated synthetic peptide derived from human RNF109/TRIM56: 651-755/755.

Purification: affinity purified by Protein A

Storage: 0.01M TBS (pH7.4) with 1% BSA, 0.02% Proclin300 and 50% Glycerol.


Shipped at 4°C. Store at -20°C for one year. Avoid repeated freeze/thaw cycles.

Background: E3 ubiquitin-protein ligase that mediates 'Lys-63'-linked polyubiquitination of


TMEM173/STING, thereby playing a key role in innate immunity. TMEM173/STING 'Lys-63'-linked ubiquitination activates the production of type I interferon IFN-beta following

detection of pathogen- and host-derived double-stranded DNA.

## **VALIDATION IMAGES**

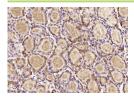

Paraformaldehyde-fixed, paraffin embedded (rat thyroid gland); Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15min; Block endogenous peroxidase by 3% hydrogen peroxide for 20 minutes; Blocking buffer (normal goat serum) at 37°C for 30min; Antibody incubation with (RNF109 TRIM56) Polyclonal Antibody, Unconjugated (bs-17116R) at 1:200 overnight at 4°C, followed by operating according to SP Kit(Rabbit) (sp-0023) instructions and DAB staining.

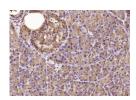

Paraformaldehyde-fixed, paraffin embedded (rat pancreas); Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15min; Block endogenous peroxidase by 3% hydrogen peroxide for 20 minutes; Blocking buffer (normal goat serum) at 37°C for 30min; Antibody incubation with (RNF109 TRIM56) Polyclonal Antibody, Unconjugated (bs-17116R) at 1:200 overnight at 4°C, followed by operating according to SP Kit(Rabbit) (sp-0023) instructions and DAB staining.