bs-5145R

[Primary Antibody]

www.bioss.com.cn sales@bioss.com.cn techsupport@bioss.com.cn 400-901-9800

proBNP Rabbit pAb

DATASHEET -

Host: Rabbit Isotype: IgG

Clonality: Polyclonal

Target: proBNP

Immunogen: KLH conjugated synthetic peptide derived from mouse proBNP:

21-121/121

Purification: affinity purified by Protein A

Concentration: 1mg/ml

Storage: 0.01M TBS (pH7.4) with 1% BSA, 0.02% Proclin300 and 50%

Shipped at 4°C. Store at -20°C for one year. Avoid repeated

freeze/thaw cycles.

Background: Brain natriuretic peptide (BNP), also known as B-type natriuretic peptide, is a hormone secreted by cardiomyocytes in the heart ventricles in response to stretching caused by increased

ventricular blood volume.

The 32-amino acid polypeptide BNP is secreted attached to a 76-amino acid N-terminal fragment in the prohormone called NTproBNP (BNPT), which is biologically inactive. Once released, BNP binds to and activates the atrial natriuretic factor receptor NPRA, and to a lesser extent NPRB, in a fashion similar to atrial natriuretic peptide (ANP) but with 10-fold lower affinity. The biological halflife of BNP, however, is twice as long as that of ANP, and that of NTproBNP is even longer, making these peptides better targets than ANP for diagnostic blood testing.

Applications: IHC-P (1:100-500)

IHC-F (1:100-500) **IF** (1:100-500)


Reactivity: Mouse, Rat

Predicted


5/11 kDa MW.:

Subcellular Location: Secreted

VALIDATION IMAGES

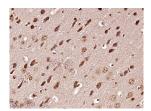

Paraformaldehyde-fixed, paraffin embedded (mouse brain tissue); Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15min; Block endogenous peroxidase by 3% hydrogen peroxide for 20 minutes; Blocking buffer (normal goat serum) at 37°C for 30min; Antibody incubation with (proBNP) Polyclonal Antibody, Unconjugated (bs-5145R) at 1:400 overnight at 4°C, followed by operating according to SP Kit(Rabbit) (sp-0023) instructions and DAB staining.

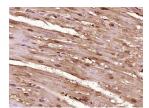

Paraformaldehyde-fixed, paraffin embedded (rat heart tissue); Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15min; Block endogenous peroxidase by 3% hydrogen peroxide for 20 minutes; Blocking buffer (normal goat serum) at 37°C for 30min; Antibody incubation with (proBNP) Polyclonal Antibody, Unconjugated (bs-5145R) at 1:400 overnight at 4°C, followed by operating according to SP Kit(Rabbit) (sp-0023) instructions and DAB staining.