

www.bioss.com.cn sales@bioss.com.cn techsupport@bioss.com.cn 400-901-9800

## Cathepsin H Rabbit pAb

Catalog Number: bs-1670R

Target Protein: Cathepsin H

Concentration: 1mg/ml

Form: Liquid

Host: Rabbit

Clonality: Polyclonal

Isotype: IgG

Applications: IHC-P (1:100-500), IHC-F (1:100-500), IF (1:100-500)

Reactivity: Human, Mouse (predicted:Rat, Rabbit, Pig, Cow, Dog, Horse)

Predicted MW: 8/35 kDa

Entrez Gene: 1512

Swiss Prot: P09668

Source: KLH conjugated synthetic peptide derived from human Cathepsin H Activation peptide:

71-170/335.

Purification: affinity purified by Protein A

Storage: 0.01M TBS (pH7.4) with 1% BSA, 0.02% Proclin300 and 50% Glycerol.


Shipped at 4°C. Store at -20°C for one year. Avoid repeated freeze/thaw cycles.

Background: The cathepsin family of proteolytic enzymes contains several diverse classes of proteases.


The cysteine protease class comprises cathepsins B, L, H, K, S, and O. The aspartyl protease class is composed of cathepsins D and E. Cathepsin G is in the serine protease class. Most cathepsins are lysosomal and each is involved in various cellular events such as peptide biosynthesis and protein degradation. The propeptide region of cathepsin H is highly homologous to other cathepsins. This region is thought to play a role in the processing of cysteine proteases. Cathepsin H is thought to be involved in several pathological conditions,

including neuro degenerative disorders and cancer.

## **VALIDATION IMAGES**

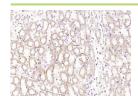

Paraformaldehyde-fixed, paraffin embedded (mouse kidney); Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15min; Block endogenous peroxidase by 3% hydrogen peroxide for 20 minutes; Blocking buffer (normal goat serum) at 37°C for 30min; Antibody incubation with (Cathepsin H) Polyclonal Antibody, Unconjugated (bs-1670R) at 1:200 overnight at 4°C, followed by operating according to SP Kit(Rabbit) (sp-0023) instructions and DAB staining.

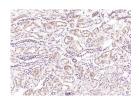

Paraformaldehyde-fixed, paraffin embedded (Human kidney); Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15min; Block endogenous peroxidase by 3% hydrogen peroxide for 20 minutes; Blocking buffer (normal goat serum) at 37°C for 30min; Antibody incubation with (Cathepsin H) Polyclonal Antibody, Unconjugated (bs-1670R) at 1:200 overnight at 4°C, followed by operating according to SP Kit(Rabbit) (sp-0023) instructions and DAB staining.