bsm-41212M

## [ Primary Antibody ]

Isotype: IgG

CloneNo.: 3A9

SWISS: P19429

## Cardiac Troponin I/TNIC Mouse mAb



www.bioss.com.cn sales@bioss.com.cn techsupport@bioss.com.cn 400-901-9800

Applications: IHC-P (1:100-500) IHC-F (1:100-500) IF (1:100-500)

Reactivity: Human, Mouse, Rat

Predicted 23 kDa MW.:

Subcellular Location: Cytoplasm

Concentration: 1mg/ml

- DATASHEET -

Host: Mouse

**Clonality:** Monoclonal


Target: Cardiac Troponin I/TNIC Purification: affinity purified by Protein A

GenelD: 7137


Storage: Size : 50ul/100ul/200ul 0.01M TBS (pH7.4) with 1% BSA, 0.02% Proclin300 and 50% Glycerol. Size : 200ug (PBS only) 0.01M PBS Shipped at 4°C. Store at -20°C for one year. Avoid repeated freeze/thaw cycles.

Background: Troponin I (TnI), along with troponin T (TnT) and troponin C (TnC), is one of 3 subunits that form the troponin complex of the thin filaments of striated muscle. Tnl is the inhibitory subunit; blocking actin-myosin interactions and thereby mediating striated muscle relaxation. The TnI subfamily contains three genes: TnI-skeletalfast-twitch, TnI-skeletal-slow-twitch, and TnI-cardiac. This gene encodes the TnI-cardiac protein and is exclusively expressed in cardiac muscle tissues. Mutations in this gene cause familial hypertrophic cardiomyopathy type 7 (CMH7) and familial restrictive cardiomyopathy (RCM). [provided by RefSeq, Jul 2008]

## VALIDATION IMAGES

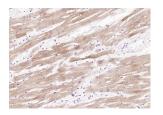

Paraformaldehyde-fixed, paraffin embedded (human heart); Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15min; Block endogenous peroxidase by 3% hydrogen peroxide for 20 minutes; Blocking buffer (normal goat serum) at 37°C for 30min; Incubation with (Cardiac Troponin I/TNIC) Monoclonal Antibody. Unconjugated (bsm-41212M) at 1:200 overnight at 4°C, followed by operating according to SP Kit(Mouse)(sp-0024) instructionsand DAB staining.

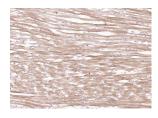

Paraformaldehyde-fixed, paraffin embedded (rat heart); Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15min; Block endogenous peroxidase by 3% hydrogen peroxide for 20 minutes; Blocking buffer (normal goat serum) at 37°C for 30min; Incubation with (Cardiac Troponin I/TNIC) Monoclonal Antibody, Unconjugated (bsm-41212M) at 1:200 overnight at 4°C, followed by operating according to SP Kit(Mouse)(sp-0024) instructionsand DAB staining.