bsm-41212M

[Primary Antibody]

Bioss ANTIBODIES

www.bioss.com.cn sales@bioss.com.cn techsupport@bioss.com.cn 400-901-9800

Cardiac Troponin I/TNIC Mouse mAb

- DATASHEET -

Host: Mouse Isotype: IgG
Clonality: Monoclonal CloneNo.: 3A9
GeneID: 7137 SWISS: P19429

Target: Cardiac Troponin I/TNIC **Purification:** affinity purified by Protein A

Concentration: 1mg/ml

Storage: Size: 50ul/100ul/200ul

0.01M TBS (pH7.4) with 1% BSA, 0.02% Proclin300 and 50%

Glycerol.

Size: 200ug (PBS only)

0.01M PBS

Shipped at 4°C. Store at -20°C for one year. Avoid repeated

freeze/thaw cycles.

Background: Troponin I (TnI), along with troponin T (TnT) and troponin C (TnC),

is one of 3 subunits that form the troponin complex of the thin filaments of striated muscle. Tnl is the inhibitory subunit; blocking actin-myosin interactions and thereby mediating striated muscle relaxation. The Tnl subfamily contains three genes: Tnl-skeletal-fast-twitch, Tnl-skeletal-slow-twitch, and Tnl-cardiac. This gene encodes the Tnl-cardiac protein and is exclusively expressed in cardiac muscle tissues. Mutations in this gene cause familial hypertrophic cardiomyopathy type 7 (CMH7) and familial restrictive cardiomyopathy (RCM). [provided by RefSeq, Jul 2008]

Applications: IHC-P (1:100-500)


IHC-F (1:100-500) IF (1:100-500)

Reactivity: Human, Mouse, Rat

Predicted MW.: ^{23 kDa}

Subcellular Cytoplasm

VALIDATION IMAGES

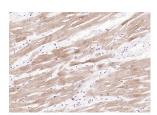

Paraformaldehyde-fixed, paraffin embedded (human heart); Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15min; Block endogenous peroxidase by 3% hydrogen peroxide for 20 minutes; Blocking buffer (normal goat serum) at 37°C for 30min; Incubation with (Cardiac Troponin I/TNIC) Monoclonal Antibody, Unconjugated (bsm-41212M) at 1:200 overnight at 4°C, followed by operating according to SP Kit(Mouse)(sp-0024) instructions and DAB staining.



Paraformaldehyde-fixed, paraffin embedded (rat heart); Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15min; Block endogenous peroxidase by 3% hydrogen peroxide for 20 minutes; Blocking buffer (normal goat serum) at 37°C for 30min; Incubation with (Cardiac Troponin I/TNIC) Monoclonal Antibody, Unconjugated (bsm-41212M) at 1:200 overnight at 4°C, followed by operating according to SP Kit(Mouse)(sp-0024) instructions and DAB staining.