



Cardiac Troponin I/TNIC Mouse mAb

Catalog Number: bsm-41212M

Target Protein: Cardiac Troponin I/TNIC

Concentration: 1mg/ml

Form: Size: 50ul/100ul/200ul

Liquid

Size: 200ug (PBS only)

Lyophilized

Note: Centrifuge tubes before opening. Reconstitute the lyophilized product in distilled

water. Optimal concentration should be determined by the end user.

Host: Mouse

Clonality: Monoclonal

Isotype: IgG

Applications: IHC-P (1:100-500), IHC-F (1:100-500), IF (1:100-500)

Reactivity: Human, Mouse, Rat

Predicted MW: 23 kDa

Subcellular Cytoplasm

Locations:

Entrez Gene: 7137 Swiss Prot: P19429

3W133 1 TOt. 1 19429

Purification: affinity purified by Protein A

Storage: Size:50ul/100ul/200ul

0.01 M TBS (pH7.4) with 1% BSA, 0.02% Proclin300 and 50% Glycerol.

Size: 200ug (PBS only)

0.01M PBS


Shipped at 4°C. Store at -20°C for one year. Avoid repeated freeze/thaw cycles.

Background: Troponin I (TnI), along with troponin T (TnT) and troponin C (TnC), is one of 3 subunits that

form the troponin complex of the thin filaments of striated muscle. TnI is the inhibitory subunit; blocking actin-myosin interactions and thereby mediating striated muscle relaxation. The TnI subfamily contains three genes: TnI-skeletal-fast-twitch, TnI-skeletal-slow-twitch, and TnI-cardiac. This gene encodes the TnI-cardiac protein and is exclusively expressed in cardiac muscle tissues. Mutations in this gene cause familial hypertrophic cardiomyopathy type 7 (CMH7) and familial restrictive cardiomyopathy (RCM). [provided by

RefSeq, Jul 2008]


VALIDATION IMAGES

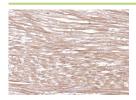

Paraformaldehyde-fixed, paraffin embedded (rat heart); Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15min; Block endogenous peroxidase by 3% hydrogen peroxide for 20 minutes; Blocking buffer (normal goat serum) at 37°C for 30min; Incubation with (Cardiac Troponin I/TNIC) Monoclonal Antibody, Unconjugated (bsm-41212M) at 1:200 overnight at 4°C, followed by operating according to SP Kit(Mouse)(sp-0024) instructions and DAB staining.



Paraformaldehyde-fixed, paraffin embedded (human heart); Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15min; Block endogenous peroxidase by 3% hydrogen peroxide for 20 minutes; Blocking buffer (normal goat serum) at 37°C for 30min; Incubation with (Cardiac Troponin I/TNIC) Monoclonal Antibody, Unconjugated (bsm-41212M) at 1:200 overnight at 4°C, followed by operating according to SP Kit(Mouse)(sp-0024) instructions and DAB staining.

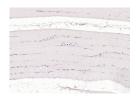


Negative sample: Paraformaldehyde-fixed, paraffin embedded (mouse skeletal muscle); Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15min; Paraformaldehyde-fixed, paraffin embedded (mouse skeletal muscle); Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15min; Block endogenous peroxidase by 3% hydrogen peroxide for 20 minutes; Blocking buffer (normal goat serum) at 37°C for 30min; Incubation with (Cardiac Troponin I/TNIC) Monoclonal Antibody, Unconjugated (bsm-41212M) at 1:200 overnight at 4°C, followed by operating according to SP Kit(Mouse)(sp-0024) instructions and DAB staining.



Negative sample: Paraformaldehyde-fixed, paraffin embedded (rat skeletal muscle);Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15min;Paraformaldehyde-fixed, paraffin embedded (rat skeletal muscle); Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15min; Block endogenous peroxidase by 3% hydrogen peroxide for 20 minutes; Blocking buffer (normal goat serum) at 37°C for 30min; Incubation with (Cardiac Troponin I/TNIC) Monoclonal Antibody, Unconjugated (bsm-41212M) at 1:200 overnight at 4°C, followed by operating according to SP Kit(Mouse)(sp-0024) instructions and DAB staining.



Negative sample: Paraformaldehyde-fixed, paraffin embedded (human skeletal muscle); Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15min; Paraformaldehyde-fixed, paraffin embedded (human skeletal muscle); Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15min; Block endogenous peroxidase by 3% hydrogen peroxide for 20 minutes; Blocking buffer (normal goat serum) at 37°C for 30min; Incubation with (Cardiac Troponin I/TNIC) Monoclonal Antibody, Unconjugated (bsm-41212M) at 1:200 overnight at 4°C, followed by operating according to SP Kit(Mouse)(sp-0024) instructions and DAB staining.

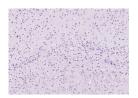

Negative sample: Paraformaldehyde-fixed, paraffin embedded (mouse brain); Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15min; Paraformaldehyde-fixed, paraffin embedded (mouse brain); Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15min; Block endogenous peroxidase by 3% hydrogen peroxide for 20 minutes; Blocking buffer (normal goat serum) at 37°C for 30min; Incubation with (Cardiac Troponin I/TNIC) Monoclonal Antibody, Unconjugated (bsm-41212M) at 1:200 overnight at 4°C, followed by operating according to SP Kit(Mouse)(sp-0024) instructions and DAB staining.