

www.bioss.com.cn sales@bioss.com.cn techsupport@bioss.com.cn 400-901-9800

DPP1 Rabbit pAb

Catalog Number: bs-2672R

Target Protein: DPP1
Concentration: 1mg/ml

Form: Liquid

Host: Rabbit

Clonality: Polyclonal

Isotype: IgG

Applications: IHC-P (1:100-500), IHC-F (1:100-500), IF (1:100-500)

Reactivity: Rat (predicted: Human, Mouse, Rabbit, Pig, Cow, Chicken, Dog, Horse)

Predicted MW: 8/52 kDa

Entrez Gene: 1075

Swiss Prot: P53634

Source: KLH conjugated synthetic peptide derived from human cathepsin C light chain: 365-463/463.

Purification: affinity purified by Protein A

Storage: 0.01M TBS (pH7.4) with 1% BSA, 0.02% Proclin300 and 50% Glycerol.

Shipped at 4°C. Store at -20°C for one year. Avoid repeated freeze/thaw cycles.

Background: The protein encoded by this gene, a member of the peptidase C1 family, is a lysosomal

cysteine proteinase that appears to be a central coordinator for activation of many serine proteinases in immune/inflammatory cells. It is composed of a dimer of disulfide-linked heavy and light chains, both produced from a single protein precursor, and a residual portion of the propeptide acts as an intramolecular chaperone for the folding and


stabilization of the mature enzyme. This enzyme requires chloride ions for activity and can

degrade glucagon. Defects in the encoded protein have been shown to be a cause of

Papillon-Lefevre syndrome, an autosomal recessive disorder characterized by palmoplantar keratosis and periodontitis. Multiple transcript variants encoding different isoforms have

been found for this gene. [provided by RefSeq].

VALIDATION IMAGES

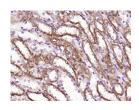

Paraformaldehyde-fixed, paraffin embedded (Rat kidney); Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15min; Block endogenous peroxidase by 3% hydrogen peroxide for 20 minutes; Blocking buffer (normal goat serum) at 37°C for 30min; Antibody incubation with (DPP1) Polyclonal Antibody, Unconjugated (bs-2672R) at 1:400 overnight at 4°C, followed by a conjugated secondary antibody (sp-0023) for 20 minutes and DAB staining.