bs-3223R

[Primary Antibody]

BIOSS

www.bioss.com.cn sales@bioss.com.cn techsupport@bioss.com.cn 400-901-9800

phospho-HSL (Ser863) Rabbit pAb

- DATASHEET -

Host: Rabbit **Isotype:** IgG

Clonality: Polyclonal

GenelD: 25330 **SWISS:** P15304

Target: HSL (Ser863)

Immunogen: KLH conjugated synthesised phosphopeptide derived from rat HSL

around the phosphorylation site of Ser863: RR(p-S)V.

Purification: affinity purified by Protein A

Concentration: 1mg/ml

Storage: 0.01M TBS (pH7.4) with 1% BSA, 0.02% Proclin300 and 50%

Glycerol.

Shipped at 4°C. Store at -20°C for one year. Avoid repeated

freeze/thaw cycles.

Background: HSL/LIPE is found in adipose tissue and heart, where it primarily

hydrolyzes stored triglycerides to free fatty acids. It is also found in steroidogenic tissues, where it principally converts cholesteryl esters to free cholesterol for steroid hormone production. There

are two named isoforms.

Applications: IHC-P (1:100-500)

IHC-F (1:100-500) **IF** (1:100-500)


Reactivity: Rat (predicted: Human,

Mouse)


Predicted MW.: 117 kDa

Subcellular Location: Cell membrane ,Cytoplasm

VALIDATION IMAGES

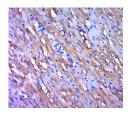

Paraformaldehyde-fixed, paraffin embedded (rat heart tissue); Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15min; Block endogenous peroxidase by 3% hydrogen peroxide for 20 minutes; Blocking buffer (normal goat serum) at 37°C for 30min; Antibody incubation with (P-HSL (Ser863)) Polyclonal Antibody, Unconjugated (bs-3223R) at 1:400 overnight at 4°C, followed by a conjugated secondary (sp-0023) for 20 minutes and DAB staining.

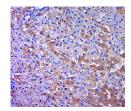

Paraformaldehyde-fixed, paraffin embedded (rat liver tissue); Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15min; Block endogenous peroxidase by 3% hydrogen peroxide for 20 minutes; Blocking buffer (normal goat serum) at 37°C for 30min; Antibody incubation with (P-HSL (Ser863)) Polyclonal Antibody, Unconjugated (bs-3223R) at 1:400 overnight at 4°C, followed by a conjugated secondary (sp-0023) for 20 minutes and DAB staining.