#### bs-3448R

### [ Primary Antibody ]

# BIOSS

www.bioss.com.cn sales@bioss.com.cn techsupport@bioss.com.cn 400-901-9800

## phospho-Tie2 (Ser1119) Rabbit pAb

- DATASHEET -

Host: Rabbit Isotype: IgG

Clonality: Polyclonal

**GeneID:** 7010 **SWISS:** Q02763

Target: Tie2 (Ser1119)

**Immunogen:** KLH conjugated Synthesised phosphopeptide derived from human

Tie2 around the phosphorylation site of Ser1119: DC(p-S)AE.

**Purification:** affinity purified by Protein A

Concentration: 1mg/ml

**Storage:** 0.01M TBS (pH7.4) with 1% BSA, 0.02% Proclin300 and 50%

Glycerol.

Shipped at 4°C. Store at -20°C for one year. Avoid repeated

freeze/thaw cycles.

Background: The TEK receptor tyrosine kinase is expressed almost exclusively in

endothelial cells in mice, rats, and humans. This receptor possesses a unique extracellular domain containing 2 immunoglobulin-like loops separated by 3 epidermal growth factor-like repeats that are connected to 3 fibronectin type III-like repeats. The ligand for the receptor is angiopoietin-1. Defects in TEK are associated with inherited venous malformations; the TEK signaling pathway appears to be critical for endothelial cell-smooth muscle cell communication in venous morphogenesis.TEK

is closely related to the TIE receptor tyrosine kinase.

Applications: IHC-P (1:100-500)

IHC-F (1:100-500) IF (1:100-500)

Reactivity: Human (predicted: Mouse,


Rat, Rabbit, Pig, Sheep, Cow, Chicken, Dog, Horse)

Predicted MW.: 124 kDa

**Subcellular** Secreted ,Extracellular **Location:** matrix ,Cell membrane

,Cytoplasm

### VALIDATION IMAGES

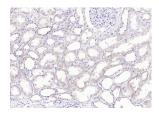

Paraformaldehyde-fixed, paraffin embedded (Human kidney); Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15min; Block endogenous peroxidase by 3% hydrogen peroxide for 20 minutes; Blocking buffer (normal goat serum) at 37°C for 30min; Antibody incubation with (Phospho-Tie2 (Ser1119))
Polyclonal Antibody, Unconjugated (bs-3448R) at 1:200 overnight at 4°C, followed by operating according to SP Kit(Rabbit) (sp-0023) instructions and DAB staining.