bs-3230R

[Primary Antibody]

BIOSS ANTIBODIES

www.bioss.com.cn sales@bioss.com.cn techsupport@bioss.com.cn 400-901-9800

phospho-IKK gamma (Ser376) Rabbit pAb

- DATASHEET -

Host: Rabbit **Isotype:** IgG

Clonality: Polyclonal

GeneID: 8517 **SWISS:** Q9Y6K9

Target: IKK gamma (Ser376)

Immunogen: KLH conjugated Synthesised phosphopeptide derived from human

IKBKG around the phosphorylation site of Ser376: YL(p-S)SP-NH2.

Purification: affinity purified by Protein A

Concentration: 1mg/ml

Storage: 0.01M TBS (pH7.4) with 1% BSA, 0.02% Proclin300 and 50%

Glycerol.

Shipped at 4°C. Store at -20°C for one year. Avoid repeated

freeze/thaw cycles.

Background: Pro inflammatory cytokines activate the transcription factor NF

kappa B by stimulating the activity of a protein kinase that phosphorylates Ikappa B, an inhibitor of NF kappa B, at sites that trigger its ubiquitination and degradation. A large, cytokine responsive Ikappa B kinase (IKK) complex has been purified and the genes encoding 2 of its subunits have been cloned. These subunits, IKK alpha and Ikk beta, are protein kinases whose function is needed for NF kappa B activation by pro inflammatory stimuli. IKK is composed of similar amounts of IKK alpha, Ikk beta, which are differentially processed forms of a third subunit, IKK gamma. IKK gamma interacts preferentially with IKK beta and is

required for the activation of the IKK complex.

Applications: IHC-P (1:100-500)

IHC-F (1:100-500) IF (1:100-500)


Reactivity: Rat (predicted: Human)

Predicted 46 kD

MW.: 46 KDa

Subcellular Cytoplasm ,Nucleus

VALIDATION IMAGES

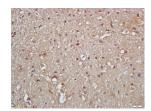

Tissue/cell: rat brain tissue; 4%
Paraformaldehyde-fixed and paraffinembedded; Antigen retrieval: citrate buffer (0.01M, pH 6.0), Boiling bathing for 15min; Block endogenous peroxidase by 3% Hydrogen peroxide for 30min; Blocking buffer (normal goat serum,C-0005) at 37°C for 20 min; Incubation: Anti-phospho-IKK gamma (Ser376) Polyclonal Antibody, Unconjugated(bs-3230R) 1:200, overnight at 4°C, followed by conjugation to the secondary antibody(SP-0023) and DAB(C-0010) staining