

bs-25423R**[Primary Antibody]****BioSS**
ANTIBODIES

www.bioss.com.cn

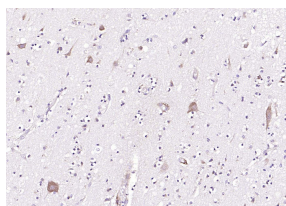
sales@bioss.com.cn

techsupport@bioss.com.cn

400-901-9800

FMN1 Rabbit pAb**— DATASHEET —**

Host: Rabbit	Isotype: IgG	Applications: IHC-P (1:100-500) IHC-F (1:400-800) IF (1:100-500) Reactivity: Human (predicted: Mouse, Rat, Rabbit, Cow, Chicken, Dog, Horse) Predicted MW.: 157 kDa Subcellular Location: Cell membrane, Cytoplasm, Nucleus
Clonality: Polyclonal		
GeneID: 342184	SWISS: Q68DA7	
Target: FMN1		
Immunogen: KLH conjugated synthetic peptide derived from human FMN1: 1141-1240/1419.		
Purification: affinity purified by Protein A		
Concentration: 1mg/ml		
Storage: 0.01M TBS (pH7.4) with 1% BSA, 0.02% Proclin300 and 50% Glycerol. Shipped at 4°C. Store at -20°C for one year. Avoid repeated freeze/thaw cycles.		
Background: The temporal genetic hierarchy influencing normal limb development can deregulate and mediate mammalian developmental syndromes. In mice, the limb deformity (ld) locus influences normal limb development and gives rise to alternative mRNAs that can translate into a family of proteins known as formins. Formins play a crucial role in cytoskeletal reorganization by influencing Actin filament assembly. Formins co-localize with the actin cytoskeleton and can translocate into the cell cytosol and into the nucleus in an HGF-dependent manner. Vertebrate nuclear formins can control polarizing activity in limb buds through establishment of a Sonic hedgehog/FGF-4 feedback loop. Deficiency mutations at the mammalian ld locus lead to profound developmental defects in limb and kidney formation. The human Formin 1 and 2 genes map to chromosome 15q13.3 and 1q43, respectively.		

— VALIDATION IMAGES —

Paraformaldehyde-fixed, paraffin embedded (human brain); Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15min; Block endogenous peroxidase by 3% hydrogen peroxide for 20 minutes; Blocking buffer (normal goat serum) at 37°C for 30min; Incubation with (FMN1) Polyclonal Antibody, Unconjugated (bs-25423R) at 1:200 overnight at 4°C, followed by operating according to SP Kit(Rabbit) (sp-0023) instructions and DAB staining.