bs-21271R

[Primary Antibody]

Bioss ANTIBODIES

www.bioss.com.cn sales@bioss.com.cn techsupport@bioss.com.cn 400-901-9800

SLC36A2 Rabbit pAb

- DATASHEET -

Host: Rabbit Isotype: IgG

Clonality: Polyclonal

GenelD: 153201 **SWISS:** Q495M3

Target: SLC36A2

Immunogen: KLH conjugated synthetic peptide derived from human SLC36A2:

261-360/483. < Extracellular >

Purification: affinity purified by Protein A

Concentration: 1mg/ml

Storage: Preservative: 0.02% Proclin300, Constituents: 1% BSA, 0.01M PBS,

pH1.4.

Shipped at 4°C. Store at -20°C for one year. Avoid repeated

freeze/thaw cycles.

Background: This gene encodes a pH-dependent proton-coupled amino acid

transporter that belongs to the amino acid auxin permease 1 protein family. The encoded protein primarily transports small amino acids such as glycine, alanine and proline. Mutations in this gene are associated with iminoglycinuria and hyperglycinuria.

[provided by RefSeq, Sep 2010]

Applications: IHC-P (1:100-500)

IHC-F (1:100-500) IF (1:100-500)


Reactivity: Mouse, Rat

(predicted: Human)


Predicted 53 kDa

Subcellular Location: Cell membrane

VALIDATION IMAGES

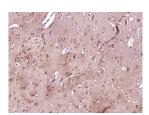

Paraformaldehyde-fixed, paraffin embedded (Rat brain); Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15min; Block endogenous peroxidase by 3% hydrogen peroxide for 20 minutes; Blocking buffer (normal goat serum) at 37°C for 30min; Antibody incubation with (SLC36A2) Polyclonal Antibody, Unconjugated (bs-21271R) at 1:400 overnight at 4°C, followed by operating according to SP Kit(Rabbit) (sp-0023) instructions and DAB staining.

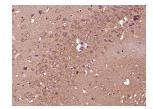

Paraformaldehyde-fixed, paraffin embedded (Mouse brain); Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15min; Block endogenous peroxidase by 3% hydrogen peroxide for 20 minutes; Blocking buffer (normal goat serum) at 37°C for 30min; Antibody incubation with (SLC36A2) Polyclonal Antibody, Unconjugated (bs-21271R) at 1:400 overnight at 4°C, followed by operating according to SP Kit(Rabbit) (sp-0023) instructions and DAB staining.