

bs-7561R**[Primary Antibody]****BioSS**
ANTIBODIES

www.bioss.com.cn


sales@bioss.com.cn

techsupport@bioss.com.cn


400-901-9800

TIRAP Rabbit pAb**— DATASHEET —**

Host: Rabbit	Isotype: IgG	Applications: IHC-P (1:100-500) IHC-F (1:100-500) IF (1:100-500) Reactivity: Human, Mouse (predicted: Rat, Rabbit) Predicted MW.: 24 kDa Subcellular Location: Cytoplasm
Clonality: Polyclonal		
GeneID: 114609	SWISS: P58753	
Target: TIRAP		
Immunogen: KLH conjugated synthetic peptide derived from human TIRAP: 81-221/221.		
Purification: affinity purified by Protein A		
Concentration: 1mg/ml		
Storage: 0.01M TBS (pH7.4) with 1% BSA, 0.02% Proclin300 and 50% Glycerol. Shipped at 4°C. Store at -20°C for one year. Avoid repeated freeze/thaw cycles.		
Background: The Toll-like receptor (TLR) family in mammals comprises a family of transmembrane proteins characterized by multiple copies of leucine rich repeats in the extracellular domain and IL-1 receptor motif in the cytoplasmic domain. Like its counterparts in Drosophila, TLRs signal through adaptor molecules. The TLR family is a phylogenetically conserved mediator of innate immunity that is essential for microbial recognition. Ten human homologs of TLRs (TLR1-10) have been described. TIRAP (TIR domain-containing adaptor protein) is an adaptor protein used by TLR4. Blocking TIRAP inhibits TLR4-mediated signaling events, including DC maturation and cytokine production.		

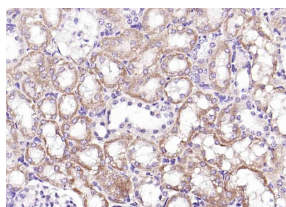— VALIDATION IMAGES —

Paraformaldehyde-fixed, paraffin embedded (mouse kidney); Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15min; Block endogenous peroxidase by 3% hydrogen peroxide for 20 minutes; Blocking buffer (normal goat serum) at 37°C for 30min; Incubation with (TIRAP) Polyclonal Antibody, Unconjugated (bs-7561R) at 1:200 overnight at 4°C, followed by operating according to SP Kit(Rabbit) (sp-0023) instructions and DAB staining.

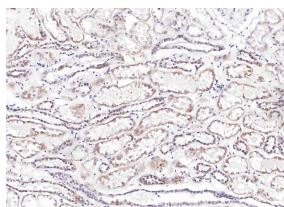

Paraformaldehyde-fixed, paraffin embedded (Human kidney); Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15min; Block endogenous peroxidase by 3% hydrogen peroxide for 20 minutes; Blocking buffer (normal goat serum) at 37°C for 30min; Incubation with (TIRAP) Polyclonal Antibody, Unconjugated (bs-7561R) at 1:200 overnight at 4°C, followed by operating according to SP Kit(Rabbit) (sp-0023) instructions and DAB staining.