

bs-13077R**[Primary Antibody]****BioSS**
ANTIBODIES

www.bioss.com.cn

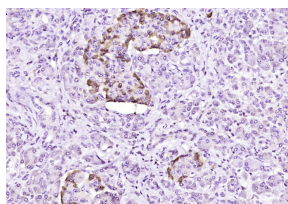
sales@bioss.com.cn

techsupport@bioss.com.cn

400-901-9800

ENT1 Rabbit pAb**— DATASHEET —**

Host: Rabbit	Isotype: IgG	Applications: IHC-P (1:100-500) IHC-F (1:100-500) IF (1:100-500)
Clonality: Polyclonal		
GeneID: 2030	SWISS: Q99808	
Target: ENT1		Reactivity: Human (predicted: Mouse, Rat, Sheep, Cow, Horse, Cynomolgus Monkey, Orangutan)
Immunogen: KLH conjugated synthetic peptide derived from human ENT1/SLC29A1: 201-300/456.		Predicted MW.: 50 kDa
Purification: affinity purified by Protein A		Subcellular Location: Cell membrane
Concentration: 1mg/ml		
Storage: 0.01M TBS (pH7.4) with 1% BSA, 0.02% Proclin300 and 50% Glycerol. Shipped at 4°C. Store at -20°C for one year. Avoid repeated freeze/thaw cycles.		
Background: Equilibrative nucleoside transporters (ENTs) regulate many physiological processes and are widely distributed in mammals, plants, yeasts, insects, nematodes and protozoans. They enable facilitated diffusion of hydrophilic nucleosides, such as adenosine and nucleoside analogs, across cell membranes. ENTs are required for uptake of antiviral and anticancer nucleoside drugs and influence a variety of physiological processes, such as neurotransmission and platelet aggregation, by regulating the amount of adenosine available to cell surface receptors. Equilibrative nucleoside transporter 1 (ENT1), also designated solute carrier family 29 (nucleoside transporters), member 1, belongs to the SLC29A transporter family and is a mammalian ENT isoform. ENT1, along with ENT3, mediates the majority of influx and efflux of nucleosides across the plasma membrane.		

— VALIDATION IMAGES —

Paraformaldehyde-fixed, paraffin embedded (human Pancreatic cancer); Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15min; Block endogenous peroxidase by 3% hydrogen peroxide for 20 minutes; Blocking buffer (normal goat serum) at 37°C for 30min; Antibody incubation with (ENT1) Polyclonal Antibody, Unconjugated (bs-13077R) at 1:200 overnight at 4°C, followed by operating according to SP Kit(Rabbit) (sp-0023) instructions and DAB staining.