

www.bioss.com.cn sales@bioss.com.cn techsupport@bioss.com.cn 400-901-9800

## **BMP7** Rabbit pAb

Catalog Number: bs-10506R

Target Protein: BMP7
Concentration: 1mg/ml

Form: Liquid Host: Rabbit

Clonality: Polyclonal

Isotype: IgG

Applications: IHC-P (1:100-500), IHC-F (1:100-500), IF (1:100-500)

Reactivity: Human (predicted:Mouse, Rat)

Predicted MW: 47 kDa Entrez Gene: 655 Swiss Prot: P18075

Purification: affinity purified by Protein A

Storage: 0.01M TBS (pH7.4) with 1% BSA, 0.02% Proclin300 and 50% Glycerol.

Shipped at 4°C. Store at -20°C for one year. Avoid repeated freeze/thaw cycles.


Background: The bone morphogenetic proteins (BMPs) are a family of secreted signaling molecules that

can induce ectopic bone growth. Many BMPs are part of the transforming growth factor-beta


(TGFB) superfamily. BMPs were originally identified by an ability of demineralized bone extract to induce endochondral osteogenesis in vivo in an extraskeletal site. Based on its expression early in embryogenesis, the BMP encoded by this gene has a proposed role in

early development and possible bone inductive activity. [provided by RefSeq].

## **VALIDATION IMAGES**

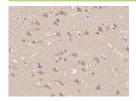

Paraformaldehyde-fixed, paraffin embedded (human brain glioma); Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15min; Block endogenous peroxidase by 3% hydrogen peroxide for 20 minutes; Blocking buffer (normal goat serum) at 37°C for 30min; Antibody incubation with (BMP7) Polyclonal Antibody, Unconjugated (bs-10506R) at 1:400 overnight at 4°C, followed by operating according to SP Kit(Rabbit) (sp-0023) instructions and DAB staining.

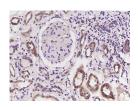

Paraformaldehyde-fixed, paraffin embedded (Human kidney tissue); Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15min; Block endogenous peroxidase by 3% hydrogen peroxide for 20 minutes; Blocking buffer (normal goat serum) at 37°C for 30min; Antibody incubation with (BMP7) Polyclonal Antibody, Unconjugated (bs-10506R) at 1:400 overnight at 4°C, followed by operating according to SP Kit(Rabbit) (sp-0023) instructions and DAB staining.