bs-16102R

[Primary Antibody]

www.bioss.com.cn sales@bioss.com.cn techsupport@bioss.com.cn 400-901-9800

AAGAB Rabbit pAb

DATASHEET -

Isotype: IgG Host: Rabbit

Clonality: Polyclonal

GenelD: 79719 SWISS: Q6PD74

Target: AAGAB

Immunogen: KLH conjugated synthetic peptide derived from human AAGAB:

221-315/315.

Purification: affinity purified by Protein A

Concentration: 1mg/ml

Storage: 0.01M TBS (pH7.4) with 1% BSA, 0.02% Proclin300 and 50%

Glycerol.

Shipped at 4°C. Store at -20°C for one year. Avoid repeated

freeze/thaw cycles.

Background: Encoding more than 700 genes, chromosome 15 is made up of approximately 106 million base pairs and is about 3% of the human genome. Angelman and Prader-Willi syndromes are associated with loss of function or deletion of genes in the 15q11q13 region. In the case of Angelman syndrome, this loss is due to inactivity of the maternal 15q11-q13 encoded UBE3A gene in the brain by either chromosomal deletion or mutation. In cases of Prader-Willi syndrome, there is a partial or complete deletion of this region from the paternal copy of chromosome 15. Tay-Sachs disease is a lethal disorder associated with mutations of the HEXA gene, which is encoded by chromosome 15. Marfan syndrome is associated with chromosome 15 through the FBN1 gene. The FLJ11506 gene product has been provisionally designated FLJ11506 pending further characterization.

Applications: IHC-P (1:100-500)

IHC-F (1:100-500) **IF** (1:100-500)

Reactivity: Rat (predicted: Human,


Mouse, Rabbit, Sheep, Cow,

Dog, Horse)

Predicted 35 kDa

Subcellular Cytoplasm

VALIDATION IMAGES

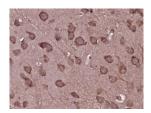

Paraformaldehyde-fixed, paraffin embedded (Rat brain); Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15min; Block endogenous peroxidase by 3% hydrogen peroxide for 20 minutes; Blocking buffer (normal goat serum) at 37°C for 30min; Antibody incubation with (AAGAB) Polyclonal Antibody. Unconjugated (bs-16102R) at 1:400 overnight at 4°C, followed by operating according to SP Kit(Rabbit) (sp-0023) instructions and DAB staining.