

bs-22754R**[Primary Antibody]****Bioss**
ANTIBODIES

www.bioss.com.cn

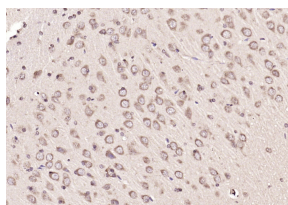
sales@bioss.com.cn

techsupport@bioss.com.cn

400-901-9800

Peroxiredoxin-5 Rabbit pAb**— DATASHEET —**

| | | |
|---|----------------------|--|
| Host: Rabbit | Isotype: IgG | Applications: IHC-P (1:100-500) |
| Clonality: Polyclonal | | IHC-F (1:100-500) |
| GeneID: 25824 | SWISS: P30044 | IF (1:100-500) |
| Target: Peroxiredoxin-5 | | |
| Immunogen: KLH conjugated synthetic peptide derived from human Peroxiredoxin-5: 51-150/214 . | | |
| Purification: affinity purified by Protein A | | Reactivity: Rat (predicted: Human, Mouse, Pig, Cow, Dog, Horse) |
| Concentration: 1mg/ml | | |
| Storage: 0.01M TBS (pH7.4) with 1% BSA, 0.02% Proclin300 and 50% Glycerol. Shipped at 4°C. Store at -20°C for one year. Avoid repeated freeze/thaw cycles. | | Predicted MW.: 22 kDa |
| Background: This gene encodes a member of the peroxiredoxin family of antioxidant enzymes, which reduce hydrogen peroxide and alkyl hydroperoxides. The encoded protein may play an antioxidant protective role in different tissues under normal conditions and during inflammatory processes. This protein interacts with peroxisome receptor 1. The crystal structure of this protein in its reduced form has been resolved to 1.5 angstrom resolution. This gene uses alternate in-frame translation initiation sites to generate mitochondrial or peroxisomal/cytoplasmic forms. Three transcript variants encoding distinct isoforms have been identified for this gene. [provided by RefSeq, Jul 2008]. | | Subcellular Location: Cytoplasm |

— VALIDATION IMAGES —

Paraformaldehyde-fixed, paraffin embedded (rat brain); Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15min; Block endogenous peroxidase by 3% hydrogen peroxide for 20 minutes; Blocking buffer (normal goat serum) at 37°C for 30min; Antibody incubation with (Peroxiredoxin-5) Polyclonal Antibody, Unconjugated (bs-22754R) at 1:200 overnight at 4°C, followed by operating according to SP Kit(Rabbit) (sp-0023) instructions and DAB staining.