## bs-22269R

## [ Primary Antibody ]

www.bioss.com.cn sales@bioss.com.cn techsupport@bioss.com.cn 400-901-9800

# AMPK alpha 2 Rabbit pAb

DATASHEET -

Host: Rabbit Isotype: IgG

Clonality: Polyclonal

**GenelD: 5563 SWISS:** P54646

Target: AMPK alpha 2

Immunogen: KLH conjugated synthetic peptide derived from human AMPK

alpha 2: 361-460/552.

Purification: affinity purified by Protein A

Concentration: 1mg/ml

Storage: 0.01M TBS (pH7.4) with 1% BSA, 0.02% Proclin300 and 50%

Glycerol.

Shipped at 4°C. Store at -20°C for one year. Avoid repeated

freeze/thaw cycles.

Background: The protein encoded by this gene is a catalytic subunit of the AMPactivated protein kinase (AMPK). AMPK is a heterotrimer consisting of an alpha catalytic subunit, and non-catalytic beta and gamma subunits. AMPK is an important energy-sensing enzyme that monitors cellular energy status. In response to cellular metabolic stresses, AMPK is activated, and thus phosphorylates and inactivates acetyl-CoA carboxylase (ACC) and beta-hydroxy betamethylglutaryl-CoA reductase (HMGCR), key enzymes involved in regulating de novo biosynthesis of fatty acid and cholesterol. Studies of the mouse counterpart suggest that this catalytic subunit may control whole-body insulin sensitivity and is necessary for maintaining myocardial energy homeostasis during ischemia. [provided by RefSeq, Jul 2008]

**Applications: IHC-P** (1:100-500)

IHC-F (1:100-500) **IF** (1:100-500)


Reactivity: Rat (predicted: Human, Pig,

Sheep, Cow, Dog, Horse)

Predicted 62 kDa MW.:

Subcellular Cytoplasm , Nucleus

### VALIDATION IMAGES

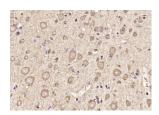

Paraformaldehyde-fixed, paraffin embedded (rat brain); Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15min; Block endogenous peroxidase by 3% hydrogen peroxide for 20 minutes; Blocking buffer (normal goat serum) at 37°C for 30min; Antibody incubation with (AMPK alpha 2) Polyclonal Antibody, Unconjugated (bs-22269R) at 1:200 overnight at 4°C, followed by operating according to SP Kit(Rabbit) (sp-0023) instructions and DAB staining.