

Phospho-POLR2A (Ser5) Ready-To-Use IHC Kit

Cat.No: IHC0626
Applications: **IHC-P**
Reactivity: Human, Mouse, Rat
Size: 50T
Assay type: Immunohistochemistry
Sample type: FFPE tissue
General Information:

Number	Component	Size	Concentration	Storage
1	PBS Buffer (powder)	2 L×2	20x	RT
2	Antigen Retrieval Buffer	20 ml	100x	2-8°C
3	Endogenous Peroxidase Blocking Buffer	3 ml	RTU	2-8°C, protect from light
4	Blocking Buffer	3 ml	RTU	2-8°C
5	Primary Antibody (Phospho-POLR2A (Ser5) Recombinant Rabbit mAb)	6 ml	RTU	2-8°C
6	Secondary Antibody (Goat Anti-Rabbit IgG H&L, HRP conjugated)	6 ml	RTU	2-8°C
7	Chromogen Component A	0.3 ml	RTU	-20°C, protect from light
8	Chromogen Component B	0.3 ml	RTU	-20°C
9	Counter Staining Reagent	5 ml	RTU	RT
10	Mounting Media	5 ml	RTU	RT
11	Control slide (Human colon, rat colon, mouse colon)	3 slides	RTU	RT
12	Datasheet	1 copy		

Storage and Stability: Please store components at the temperatures indicated on the individual tube labels. The kit is stable for 6 months from the date of receipt.

Immunohistochemistry Protocol:

1. Deparaffinization And Rehydration

Immerse slides in fresh xylene for 15 minutes and then repeat two more times using separate containers. Immerse slides sequentially in 100%, 95%, 90%, 80%, and 70% ethanol solutions for 5 minutes each. Rinse slides 3 times with distilled water for 5 minutes each.

2. Antigen Retrieval

Add 100× **Antigen Retrieval Buffer** into distilled water to prepare a 1× solution. Boil slides in 1× solution at 95°C-100°C for 15 minutes. Move the slides to 1× solution at room temperature (RT) and allow them to stand for 20 minutes. Rinse 3 times with **PBS Buffer** (dissolve the powder in 2L distilled water) for 5 minutes each.

3. **Block Endogenous Peroxidase**

Drain the liquid off the slides and then use a hydrophobic IHC pen to draw circles on the slides around tissue sections. Add 2-4 drops of **Endogenous Peroxidase Blocking Buffer** directly on slides, covering the whole tissue and block slides for 15 minutes at RT. Rinse 3 times with **PBS Buffer** for 5 minutes each.

4. **Serum Blocking**

Block with 2-4 drops of **Blocking Buffer** for 20 minutes at RT.

5. **Primary Antibody Incubation**

Drain blocking buffer from slides. Incubate slides with 2-4 drops of **Phospho-POLR2A (Ser5) Recombinant Rabbit mAb** overnight at 4°C or 1-2 hours at RT. Rinse 3 times with **PBS Buffer** for 5 minutes each.

6. **Secondary Antibody Incubation**

Incubate slides with 2-4 drops of **Goat Anti-Rabbit IgG H&L, HRP conjugated** for 1-2 hours at RT. Rinse slides 3 times with **PBS Buffer** for 5 minutes each.

7. **Signal Development**

Remove residual liquid around the tissue section. Add 50ul fresh **DAB Buffer (Chromogen Component A : Chromogen Component B : PBS Buffer=1:1:18)** to cover the tissue. Monitor the reaction under the microscope until a brown color is visible (approximate 3-5 minutes at RT). Stop reaction immediately by rinsing with distilled water. Rinse slides 3 times with distilled water for 5 minutes each.

8. **Counterstain**

Counterstain with an appropriate amount of **Counter Staining Reagent** for 3-5 minutes at RT. Rinse slides with distilled water for 5 minutes. Use 2-4 drops of **Differentiation reagent** to cover the tissue for 30 seconds. Rinse slides twice with distilled water for 5 minutes each.

9. **Dehydration Sheet**

Immerse slides sequentially in 70%, 80%, 90%, 95%, and 100% ethanol for 5 minutes each at RT. Immerse slides in 2 changes of fresh xylene, 15 minutes each. Drop some **Mounting Media** on the tissue. Mount coverslips.

Notes:

1. The positive control slide provided in the kit allows you to be sure that the experimental set-up is working properly.
2. Do not allow slides to dry at any time during this procedure.


3. Please don't replace the matching reagents in this product with other manufacturers' products.
4. As DAB is a carcinogen, please take necessary precautions.
5. PBS (reagent 1) can be stored for one week at 4°C after preparation; The antigen retrieval buffer (1× reagent 2) and the chromogenic agent (the mixture of reagents 7 and 8) should be prepared right before each assay.

Please cite this product as "IHC0626, Bioss Antibodies". Citation example: "Tissue sections using POLR2A IHC Kit (IHC0626, Bioss Antibodies) were stained for POLR2A according to the manufacturer's instructions."


Introduction:

DNA-directed RNA polymerase II subunit RPB1 (POLR1A) is a DNA-dependent RNA polymerase that catalyzes the transcription of DNA into RNA using the four ribonucleoside triphosphates as substrates. POLR1A is the largest subunit and is a catalytic component of RNA polymerase II which synthesizes mRNA precursors and many functional non-coding RNAs. It also forms the polymerase active center together with the second largest subunit. Pol II is the central component of the basal RNA polymerase II transcription machinery. It is composed of mobile elements that move relative to each other. RPB1 is part of the core element with the central large cleft, the clamp element that moves to open and close the cleft, and the jaws that are thought to grab the incoming DNA template. At the start of transcription, a single-stranded DNA template strand of the promoter is positioned within the central active site cleft of Pol II. A bridging helix emanates from RPB1 and crosses the cleft near the catalytic site and is thought to promote translocation of Pol II by acting as a ratchet that moves the RNA-DNA hybrid through the active site by switching from straight to bent conformations at each step of nucleotide addition. During transcription elongation, Pol II moves on the template as the transcript elongates. Elongation is influenced by the phosphorylation status of the C-terminal domain (CTD) of Pol II's largest subunit (RPB1), which serves as a platform for assembly of factors that regulate transcription initiation, elongation, termination, and mRNA processing. Regulation of gene expression levels depends on the balance between methylation and acetylation levels of the CTD-lysines.


Validation Data

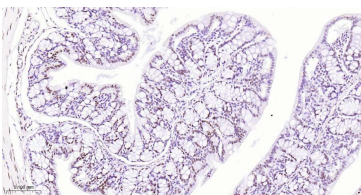

Immunohistochemical analysis of paraffin embedded Rat colon tissue slide using IHC0626 (Phospho-POLR2A (Ser5) Kit).

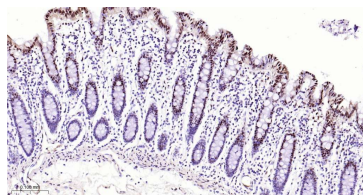

Immunohistochemical analysis of paraffin embedded Human colon tissue slide using IHC0626 (Phospho-POLR2A (Ser5) Kit).

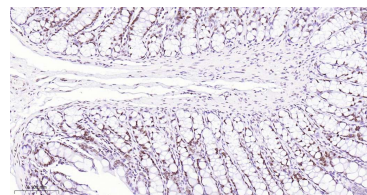

Immunohistochemical analysis of paraffin embedded Mouse colon tissue slide using IHC0626 (Phospho-POLR2A (Ser5) Kit).