bs-22112R

- DATASHEET -

[Primary Antibody]

PLK4/STK18 Rabbit pAb



www.bioss.com.cn sales@bioss.com.cn techsupport@bioss.com.cn 400-901-9800

Applications: IHC-P (1:100-500) IHC-F (1:100-500) IF (1:100-500)

Reactivity: Human (predicted: Mouse, Rat, Rabbit, Sheep, Cow)

Predicted MW.: 108 kDa

Subcellular Location: Cytoplasm ,Nucleus

Clonality: Polyclonal GenelD: 10733

Host: Rabbit

SWISS: 000444

Isotype: IgG

Target: PLK4/STK18

Immunogen: KLH conjugated synthetic peptide derived from human PLK4/STK18: 881-970/970 .


Purification: affinity purified by Protein A

Concentration: 1mg/ml

Storage: 0.01M TBS (pH7.4) with 1% BSA, 0.02% Proclin300 and 50% Glycerol. Shipped at 4°C. Store at -20°C for one year. Avoid repeated freeze/thaw cycles.

Background: Serine/threonine-protein kinase that plays a central role in centriole duplication. Able to trigger procentriole formation on the surface of the parental centriole cylinder, leading to the recruitment of centriole biogenesis proteins such as SASS6, CENPJ/CPAP, CP110, CEP135 and gamma-tubulin. When overexpressed, it is able to induce centrosome amplification through the simultaneous generation of multiple procentrioles adjoining each parental centriole during S phase. Its central role in centriole replication suggests a posible role in tumorigenesis, centrosome aberrations being frequently observed in tumors. Also involved in trophoblast differentiation by phosphorylating HAND1, leading to disrupt the interaction between HAND1 and MDFIC and activate HAND1. Phosphorylates CDC25C and CHEK2/CHK2.

- VALIDATION IMAGES

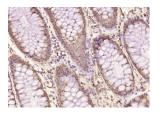

Paraformaldehyde-fixed, paraffin embedded (human colon carcinoma); Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15min; Block endogenous peroxidase by 3% hydrogen peroxide for 20 minutes; Blocking buffer (normal goat serum) at 37°C for 30min; Antibody incubation with (PLK4 STK18) Polyclonal Antibody, Unconjugated (bs-22112R) at 1:400 overnight at 4°C, followed by operating according to SP Kit(Rabbit) (sp-0023) instructionsand DAB staining.