#### bs-22140R

### [ Primary Antibody ]

# Bioss ANTIBODIES

www.bioss.com.cn sales@bioss.com.cn techsupport@bioss.com.cn 400-901-9800

**IF** (1:100-500)

## NR2F2 Rabbit pAb

- DATASHEET -

Host: Rabbit Isotype: IgG

Clonality: Polyclonal

**GenelD:** 7026 **SWISS:** P24468

Target: NR2F2

**Immunogen:** KLH conjugated synthetic peptide derived from human NR2F2:

111-210/414.

**Purification:** affinity purified by Protein A

Concentration: 1mg/ml

**Storage:** 0.01M TBS (pH7.4) with 1% BSA, 0.02% Proclin300 and 50%

Glycerol.

Shipped at 4°C. Store at -20°C for one year. Avoid repeated

freeze/thaw cycles.

**Background:** NR2F2 is involved in the regulation of the apolipoprotein A-I gene

transcription. It binds to DNA site A as a homodimer. It belongs to the nuclear hormone receptor family, NR2 subfamily. It has roles in

angiogenesis, vascular remodeling and heart development.

## 100 301 3000

**Applications: IHC-P** (1:100-500) **IHC-F** (1:100-500)


Reactivity: Human, Mouse, Rat

(predicted: Rabbit, Pig, Chicken, Dog, Horse)


Predicted MW.: 46 kDa

Subcellular Location: Nucleus

#### VALIDATION IMAGES

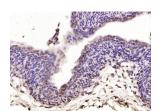

Paraformaldehyde-fixed, paraffin embedded (rat embryonic kidney); Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15min; Block endogenous peroxidase by 3% hydrogen peroxide for 20 minutes; Blocking buffer (normal goat serum) at 37°C for 30min; Antibody incubation with (NR2F2) Polyclonal Antibody, Unconjugated (bs-22140R) at 1:400 overnight at 4°C, followed by operating according to SP Kit(Rabbit) (sp-0023) instructions and DAB staining.

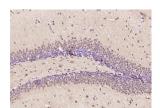

Paraformaldehyde-fixed, paraffin embedded (mouse brain); Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15min; Block endogenous peroxidase by 3% hydrogen peroxide for 20 minutes; Blocking buffer (normal goat serum) at 37°C for 30min; Antibody incubation with (NR2F2) Polyclonal Antibody, Unconjugated (bs-22140R) at 1:400 overnight at 4°C, followed by operating according to SP Kit(Rabbit) (sp-0023) instructions and DAB staining.