bs-13266R

[Primary Antibody]


Galectin 10 Rabbit pAb




www.bioss.com.cn sales@bioss.com.cn techsupport@bioss.com.cn 400-901-9800

– DATASHEET ––––		400-901-9800
Host: Rabbit	Isotype: IgG	Applications: IHC-P (1:100-500) IHC-F (1:100-500)
Clonality: Polyclonal	l	IF (1:100-500)
GenelD: 1178	SWISS: Q05315	Departivity Users and Date
Target: Galectin 1	0	Reactivity: Human, Rat
Immunogen: KLH conju 51-142/14	gated synthetic peptide derived from human Galectin 10: 2.	
Purification: affinity purified by Protein A		Predicted MW.: ^{16 kDa}
Concentration: 1mg/ml		MW.: 10 KB4
Storage: 0.01M TBS (pH7.4) with 1% BSA, 0.02% Proclin300 and 50% Glycerol. Shipped at 4°C. Store at -20°C for one year. Avoid repeated freeze/thaw cycles.		Subcellular Location: Secreted ,Cytoplasm
Background: Charcot-Leyden crystals are endogenous hexagonal bipyramidal crystals present in human tissues and secretions. Presence of Charcot-Leyden crystals correlates with the increased numbers of peripheral blood or tissue eosinophils that occurs with parasitic and allergic processes. Galectin-10, also referred to as Charcot-Leyden crystal (CLC) protein, singularly makes up these crystals. Galectin-10, a member of the galectin family of beta-Galactoside binding proteins that bind to S-type animal lectins, is expressed solely in eosinophilic and basophilic leukocytes. Galectin-10 may possess carbohydrate or IgE-binding activities, and it plays a functional role in the biology of inflammation. Expression of galectin-10 is transcriptionally induced by butyric acid.		

- VALIDATION IMAGES

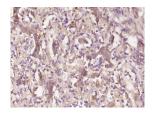

Paraformaldehyde-fixed, paraffin embedded (human liver carcinoma); Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15min; Block endogenous peroxidase by 3% hydrogen peroxide for 20 minutes; Blocking buffer (normal goat serum) at 37°C for 30min; Antibody incubation with (Galectin 10) Polyclonal Antibody, Unconjugated (bs-13266R) at 1:400 overnight at 4°C, followed by operating according to SP Kit(Rabbit) (sp-0023) instructionsand DAB staining.

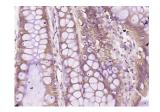

Paraformaldehyde-fixed, paraffin embedded (human colon carcinoma); Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15min; Block endogenous peroxidase by 3% hydrogen peroxide for 20 minutes; Blocking buffer (normal goat serum) at 37°C for 30min; Antibody incubation with (Galectin 10) Polyclonal Antibody, Unconjugated (bs-13266R) at 1:400 overnight at 4°C, followed by operating according to SP Kit(Rabbit) (sp-0023) instructionsand DAB staining.

Important Note: This product as supplied is intended for research use only, not for use in human, therapeutic or diagnostic applications.