bs-11219R

[Primary Antibody]

BIOSS ANTIBODIES

www.bioss.com.cn sales@bioss.com.cn techsupport@bioss.com.cn 400-901-9800

Exportin 5 Rabbit pAb

- DATASHEET -

Host: Rabbit Isotype: IgG

Clonality: Polyclonal

GenelD: 57510 SWISS: Q9HAV4

Target: Exportin 5

Immunogen: KLH conjugated synthetic peptide derived from human Exportin 5:

221-320/1204.

Purification: affinity purified by Protein A

Concentration: 1mg/ml

Storage: 0.01M TBS (pH7.4) with 1% BSA, 0.02% Proclin300 and 50%

Glycerol.

Shipped at 4°C. Store at -20°C for one year. Avoid repeated

freeze/thaw cycles.

Background: Exportin 5 (XPO5) is a member of the importin β family of proteins

and functions in tRNA export in a sequence dependent fashion. More recently, it has been shown to export pre-miRNA by a RanGTPase-driven process from the nucleus to the cytoplasm, where pre-miRNA processing occurs to produce mature miRNAs. Study of the miRNA biosynthesis pathway is essential in understanding the process of oncogenesis as global downregulation of miRNAs and the resulting alterations in the expression of tumor suppressor and oncogenic proteins is a common phenotype of cancers cells . Research studies have shown disruption of exportin 5 functions in many types of cancers including breast and lung, where pre-miRNA accumulates in the

nucleus and miRNA maturation is impaired.

Applications: IHC-P (1:100-500)

IHC-F (1:100-500) **IF** (1:100-500)


Reactivity: Mouse (predicted: Human,

Rat, Rabbit, Pig, Sheep, Cow, Chicken, Dog, Horse)

Predicted 136 kDa

Subcellular Location: Cytoplasm ,Nucleus

VALIDATION IMAGES

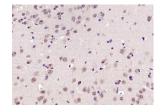

Paraformaldehyde-fixed, paraffin embedded (mouse brain); Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15min; Block endogenous peroxidase by 3% hydrogen peroxide for 20 minutes; Blocking buffer (normal goat serum) at 37°C for 30min; Antibody incubation with (Exportin 5) Polyclonal Antibody, Unconjugated (bs-11219R) at 1:200 overnight at 4°C, followed by operating according to SP Kit(Rabbit) (sp-0023) instructions and DAB staining.